

Remotely Monitor the Alcohol Level of Multiple Users From a Single Dashboard



- ✓ Access video and location verified breathalyzer results, *from any device*
- ✓ Set random, recurring or custom testing schedules for individual users
- ✓ Receive real-time email or SMS notifications for failed and missed tests



COMPANY OVERVIEW

Founded: 2001

Number of Employees: 35

Headquarters: San Francisco, CA

Other Locations: San Diego, CA, Philadelphia, PA, Miami, FL, Austin, TX

Founded in 2001, BACtrack is one of the pioneers of the breathalyzer industry. In 2004, it was the first company to obtain FDA approval to manufacture and market breathalyzers for consumer use in the United States. Since then, BACtrack has grown into a global leader in the alcohol testing and monitoring space, providing some of the most accurate and innovative breath alcohol testing devices and software on the market.


Launched in 2018, BACtrack View Enterprise was designed to fit the specific needs of treatment centers, health clinics, businesses and family law courts, making it easy and affordable to remotely monitor the sobriety of 1-1000+ users, from one centralized web portal. It currently supports thousands of users, testing multiple times per day, typically seven days a week. To date, over five million video breath alcohol tests have been recorded, processed and securely stored on the BACtrack View Enterprise platform.

ABOUT THE SERVICE



BACtrack View Enterprise (<https://enterprise.bactrack.com>) is an enterprise grade, web-based remote alcohol monitoring platform, designed to track the sobriety of users by providing a video of their breath alcohol tests to the person or organization monitoring them.

Monitored Users (MUs) are given a BACtrack Mobile Smartphone Breathalyzer to use while on the service, which connects to the BACtrack View Enterprise App on either an iPhone or Android device via Bluetooth. When it's time to take a test, the MU receives a notification on their phone, which prompts them to turn on their breathalyzer and record themselves taking the test. The BAC results and video, along with the date, time and GPS location are uploaded to the web portal in near real-time for review.

Each MU is assigned a Manager, who can request on-demand tests, set recurring or randomized testing schedules, view and approve/disapprove previous test video results, and download a summary CSV report of that user's test history. Managers have the option to receive Email or SMS notifications when an MU fails or misses a test. Additional Managers within the organization also have the option to receive these notifications.

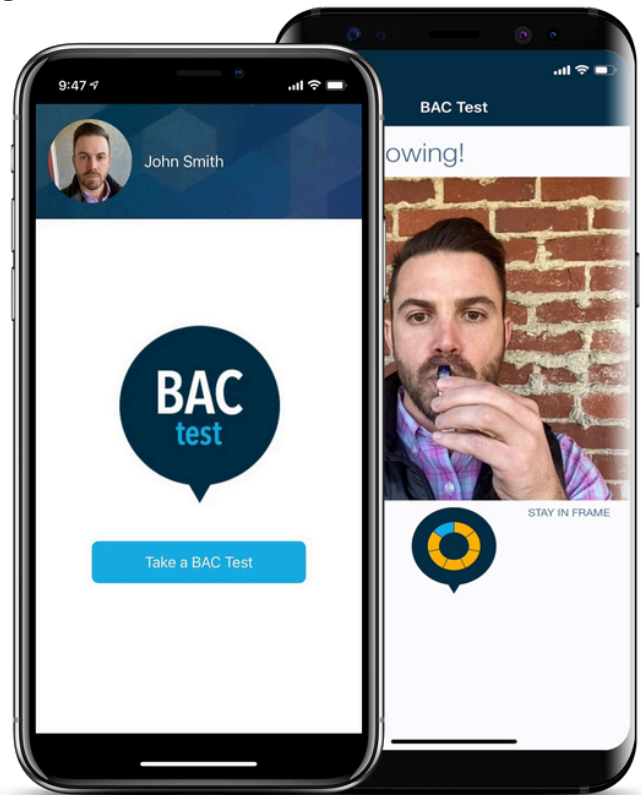

SOFTWARE SCREENSHOTS

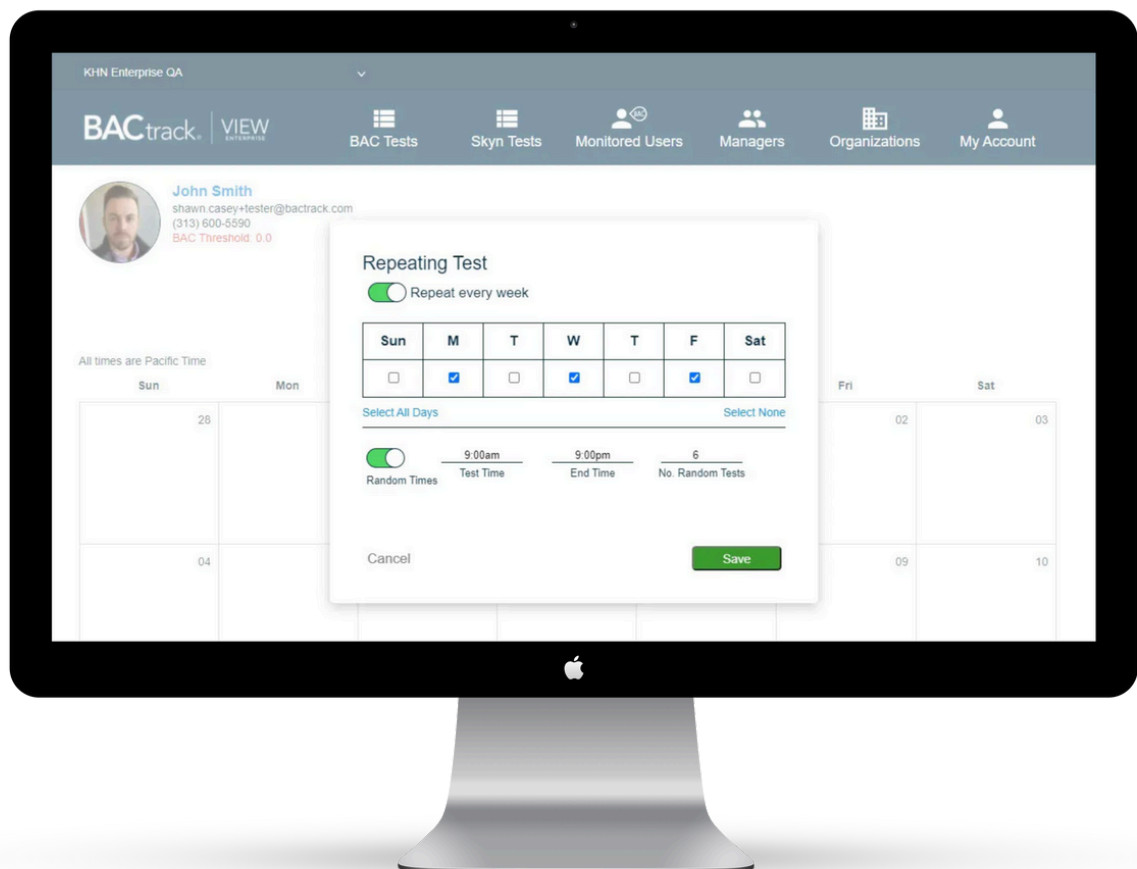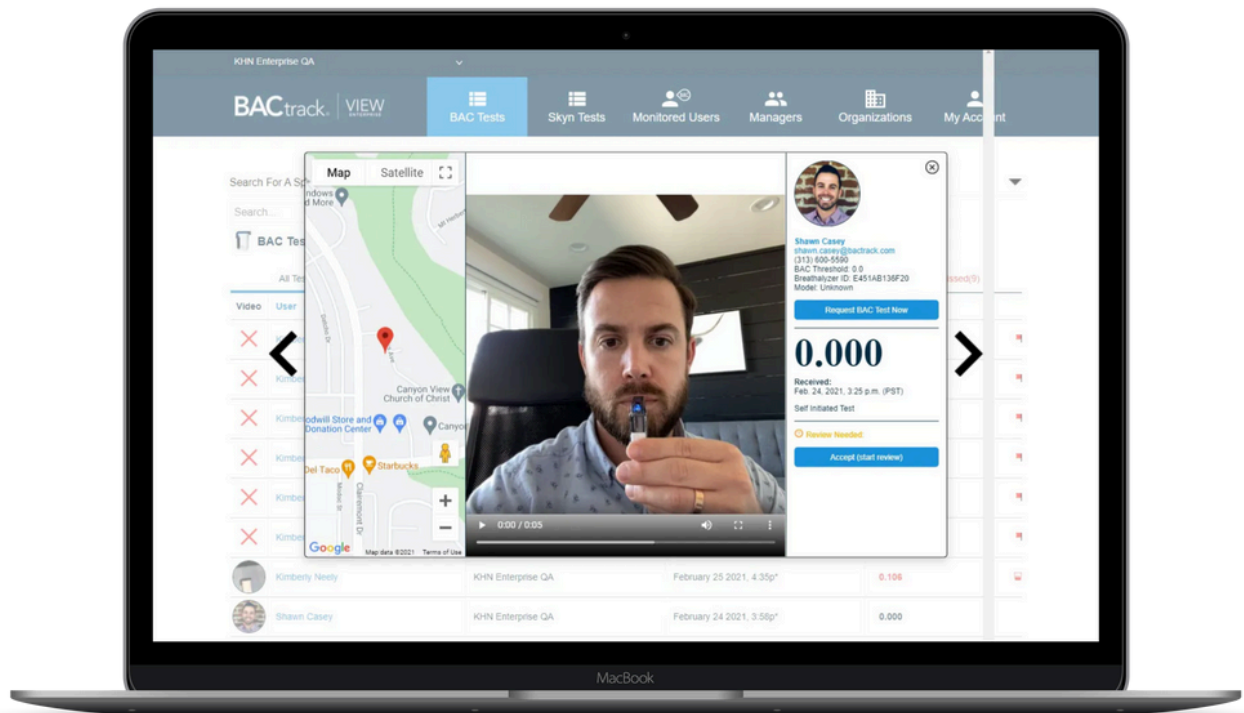

SERVICE SUMMARY

- **Web platform for Managers**

- Add new Monitored Users to the Organization
- Request on-demand BAC tests
- Set recurring or randomized testing schedules
- View and approve/disapprove video BAC results
- Download CSV testing history

- **Smartphone App for Monitored Users**

- Both iOS and Android
- Push notifications when it's time to take a test
- Records video of users during their test
- Can also take on-demand BAC tests at any point

- **High quality, discrete, fuel-cell breath alcohol tester**

- Bluetooth connectivity to App
- Charges via micro USB cable (included) - 200-300 tests on one battery charge
- 2 washable and reusable mouthpieces included, required for proper test results

CURRENT & PAST PARTNERS

- | | |
|-------------------------------|--|
| • H.A.S. of America | • Harmony Haus Sober Living |
| • The Hudson Center NY | • HeartQuest Addiction Wellness Centre |
| • Sanctum Healthcare | • Integrative Health |
| • York Children's Aid Society | • Shark Addiction Treatment |
| • Aeromedical Advisors | • Still Water Counseling |
| • Bateman Law | • Old Dominion Counseling |
| • CCSB Law | • Clear Present Recovery |
| • Suntra Modern Recovery | • Bridges Medical Clinic |
| • St. Croix County DHHS | • Alaska Pretrial Services |
| • VIP Monitoring Services | • Advanced Community Solutions |
| • NCare Health | • Track Group |

PRICING

Number of Users	Monthly Cost Per User	3 Month Pre-Paid Cost Per User	6 Month Pre-Paid Cost Per User	12 Month Pre-Paid Cost Per User
1-49	\$99.99	\$249.99	\$499.99	\$999.99
50-149	\$79.99	\$199.99	\$399.99	\$799.99
150-249	\$59.99	\$149.99	\$299.99	\$599.99
250-499	\$49.99	\$124.99	\$249.99	\$499.99
500-999	\$39.99	\$99.99	\$199.99	\$399.99
1000+	\$34.99	\$87.99	\$174.99	\$349.99